

ACHIEVE TRUE 360°

WELLNESS



THERALIGHT 360i

SAY HELLO TO THE MOST
ADVANCED FULL-BODY
PHOTOBIMODULATION
SYSTEM IN THE INDUSTRY.

Backed by years of extensive engineering and clinical research, the TheraLight 360i has evolved to become the most advanced, Versatile and effective full-body light therapy system on the market. The TheraLight 360i utilizes exclusive photobiomodulation technologies to deliver maximum tissue saturation, regardless of skin color or body type. By combining beautiful design, cutting-edge technology and precision craftsmanship, TheraLight 360i has set the standard for full-body light therapy systems.



❖ WHAT IS PHOTOBIO-MODULATION?

Photobiomodulation therapy (PBMT) is the application of red and near-infrared light to tissue that has been compromised by injury, disease or degeneration. This application increases circulation and cellular energy production while reducing inflammation and oxidative stress. PBMT increases the rate and quality of tissue repair, improves muscle performance, enhances recovery and reduces pain. PBMT is a safe and relaxing mode of recovery that has no known negative side effects.



❖ HOW DOES IT WORK?

The TheraLight 360i delivers PBMT via calibrated and adjustable penetrating wavelengths at the cellular level, which increases energy production and circulation within the body's tissues. More than 500 published double-blind, randomly controlled studies and 5,000 laboratory studies have been performed studying the effects of PBMT. The consensus shows that PBMT is a dose-dependent technology that, when properly calibrated, can deliver profound effects on a wide spectrum of conditions which include:

- Arthritis
- Tendinopathies
- Nerve pain
- Soft tissue healing
- Nerve regeneration
- Musculoskeletal pain
- Anxiety and depression
- Sleep disturbances



❖ HEAD AND FOOT CAPS

The TheraLight 360i stands as the first and only Red Light Therapy Bed to offer a TRUE 360 degree treatment. Our industry leading 45,000 diodes now adds head and foot panels allowing TheraLight 360i to offer a more direct penetration of light on the top of the head and bottom of the feet. With the capability to flip each panel up and down depending on client comfort and choice, it allows another avenue for a more advanced and customized treatment.

❖ PEMF INTEGRATION

Now offering the capability to bio-stack, the 360i is the only product within the line to integrate PEMF. The combination of PBMT and PEMF present the consumer an opportunity to target cellular repair from more than one modality at the same time

❖ APPLE IOS TECHNOLOGY

Each TheraLight 360i will now be controlled with apple ISO technology



THERALIGHT 360i

PRODUCT SPECIFICATIONS

MODEL	▶ THERALIGHT 360i
MSRP	▶ \$89,997
WARRANTY	▶ 5 YEARS
LED TYPES	▶ DAAIAS DIODE
LED RATING	▶ 130,000 HOURS
WAVELENGTHS	▶ 633nm, 810nm, 850nm, 940nm
IRRADIANCE	▶ 165 MW/CM2 +/- 20% TOLERANCE
OPERATION MODE	▶ CW AND VARIABLE
PULSED SETTINGS	▶ 0-5000HZ IN 1, 10 OR 100MH
VOLTAGE	▶ 220V / 30AMP SINGLE PHASE
INTERNAL CONTROL UNIT	▶ 7" BUILT IN TABLET
EXTERNAL CONTROL UNIT	▶ 10" LCD TOUCH SCREEN TABLET
RAPID TREATMENT TIME	▶ YES
HIGHLY REFLECTIVE PANELS	▶ YES
DIODES	▶ 45,000+
OPTIONAL PEMF	▶ YES
HEAD AND FOOT CAPS	▶ YES



CONTROL TABLET FEATURES:

- APPLE ISO APP FOR EXTERNAL TABLET
- 7" BUILT IN TOUCH SCREEN
- FULL-COLOR USER-FRIENDLY INTERFACE
- PRESET PROTOCOLS
- ABILITY TO ADJUST TIME
- INDEPENDENT CONTROL OF EACH WAVELENGTH POWER AND HZ
- ADDITIONAL PROTOCOL CUSTOMIZATION

OVERALL DIMENSIONS:

7'3"L x 4'0"W x 5'4"H
(220 cm x 122 cm x 163 cm)

TREATMENT BED LENGTH:

7'1" (211 cm)

WEIGHT:

620 lbs (281 kg)

WEIGHT CAPACITY:

500 lbs (227 kg)



THERALIGHT 360i

PEMF INTEGRATION SPECIFICATIONS

POWER LEVEL	▶	UP TO 40 GAUSS
REPEAT WAVEFORM	▶	YES
PULSED WAVEFORM	▶	YES
PRE-SET PROTOCOLS	▶	YES
ADJUSTABLE SETTING	▶	YES
MULTI - LAYERED APPROACH	▶	YES
MSRP	▶	\$10,000


Each TheraLight 360i Bed offers the ability to integrate PEMF, also known as Pulsed Electromagnetic Field.

⚡ PEMF INTEGRATION

The TheraLight 360i PEMF system utilizes both repeat waveforms and pulsed signal, strategically optimizing maximum therapeutic benefits with the body. With 4 preset protocols, Sleep Mode, Recovery Mode, Immunity Mode and Energy Mode, consumers have the ability to choose a pre-set gauss level or customize themselves up to an FDA maximum of 40Gauss

Generic PEMF devices repeatedly use the same frequency, causing cells to adapt and reduce effectiveness over time. Our PEMF dynamically adjusts frequencies, maintaining consistent therapeutic benefits without diminishing returns.

Most PEMF devices apply the same setting everywhere, our system delivers customized frequencies and intensities for each body region, ensuring precise, measurable therapeutic outcomes.

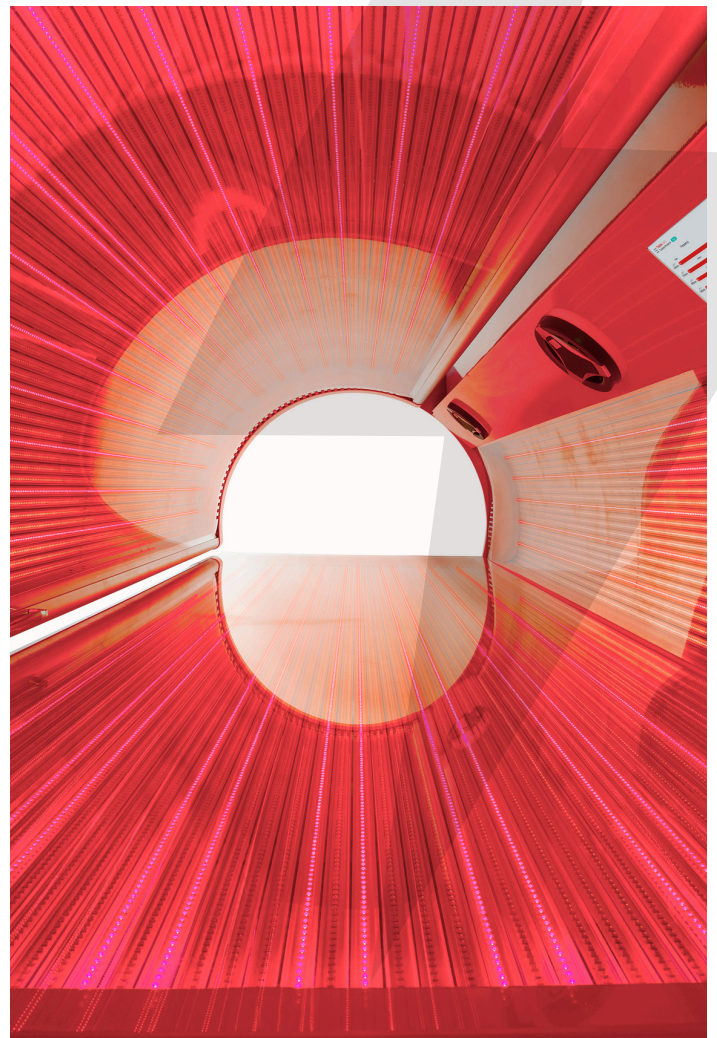



The TheraLight 360i Full Body Wellness System is registered as both a Class I and Class II device with the FDA.

Class I intended uses:

- Restoration of motion to joints
- Redevelopment of muscles
- Increased weight loss as an adjunct to a diet and exercise program

Class II intended uses:

- Relaxation of muscles and relief from muscle spasms
- Temporary relief of minor muscle and joint aches, pain and stiffness
- Temporary relief of minor pain and stiffness associated with arthritis
- Temporary increase in blood circulation



THERALIGHT LLC

TheraLight.com | 877.782.7736

175 1800 W
Suite 106 -108
Lindon UT 84042

©2025, TheraLight LLC, an FDA-Registered Emerging Leader in Photomedicine